

LINEARZ

Powered by high-intensity focused ultrasound energy, it supports both Linear Type and Dot Type with a single cartridge, enabling precise targeting of various skin depths.

- 2 **Types Mode**
Linear & Dot
- 4 **Cartridges**
Face & Body
- 7 **Depth Adjustment**
With a Single Cartridge



LINEARZ

Specifications

| | |
|--------------|---|
| Output Power | 0.1 - 3.0J |
| Spacing | Dot mode : 1.0 - 2.0 mm (0.1 mm step) Linear mode : N/A(CW) |
| Length | Dot mode : 5.0 - 25 mm (5 mm step) Linear mode : 25mm fixed |
| Cartridge | Basic (Dot/Linear) : 1.5-3.0mm (step: 0.5mm) Essential (Dot/Linear) : 1.5-3.5mm (step: 0.5mm) Core (Dot/Linear) : 3.0-6.0mm (step: 0.5mm) Contour (Dot/Linear) : 9.0mm / 11.0mm / 13.0mm |
| Dimensions | 1207(H) x 615.5(L) x 500(W) mm |
| Weight | 51 kg |



Scan or visit to learn more

Jeisys

Jeisys Medical Inc. 307, 308, 401, 410 Daeryung Techno Town 8th, Gamasanh-ro 96, Geumcheon-gu, Seoul, Korea, 08501
Website : www.jeisys.com Tel : +82 2 2603 6417 Fax : +82 2 2603 6447 Email : sales@jisys.com [H01A-C-802]

Jeisys

LINEARZ



LINEARZ

- Linear & Dot modes without changing cartridges
- Adjustable treatment depth within a single cartridge
- Versatile cartridges designed for face and body lifting



Line Lifting

HIFU lifting that penetrates deep into the fat layer



Contour Cartridge




Start customized lifting tailored to specific skin layers and area concerns.

| Basic | Essential | Core | Contour |
|------------------------------------|------------------------------------|------------------------------------|-----------------------------------|
| 4 Depths | 5 Depths | 7 Depths | 3 Depths |
| 1.5mm ~ 3.0mm 0.5mm Step / 7MHz | 1.5mm ~ 3.5mm 0.5mm Step / 7MHz | 3.0mm ~ 6.0mm 0.5mm Step / 4MHz | 9.0mm ~ 13mm 2.0mm Step / 2MHz |


Depth Adjustment

Delivers energy to specific skin layers using 4 cartridges capable of adjusting energy type and depth.

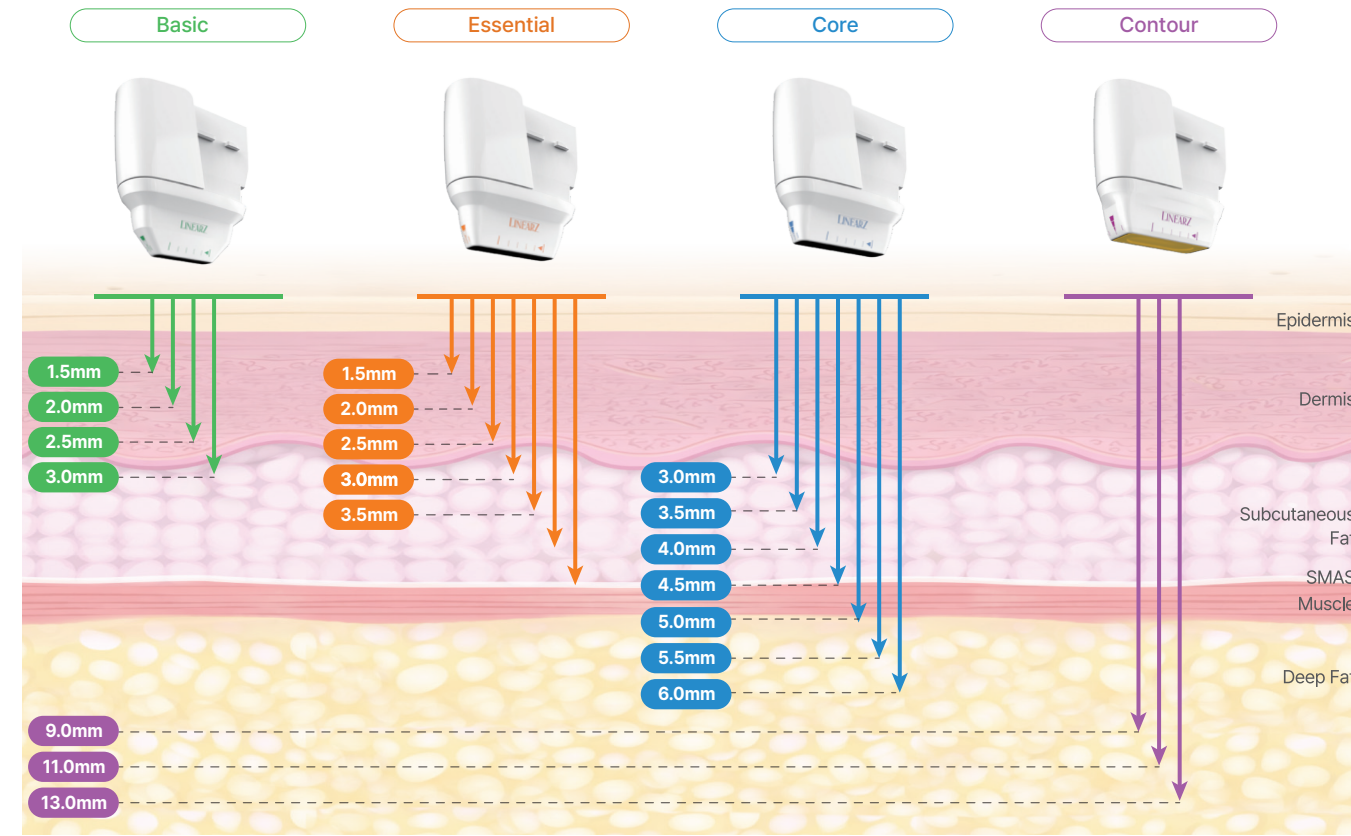

Dual-Mode Switching

Switch modes without changing cartridges.
DOT MODE <<>> LINEAR MODE

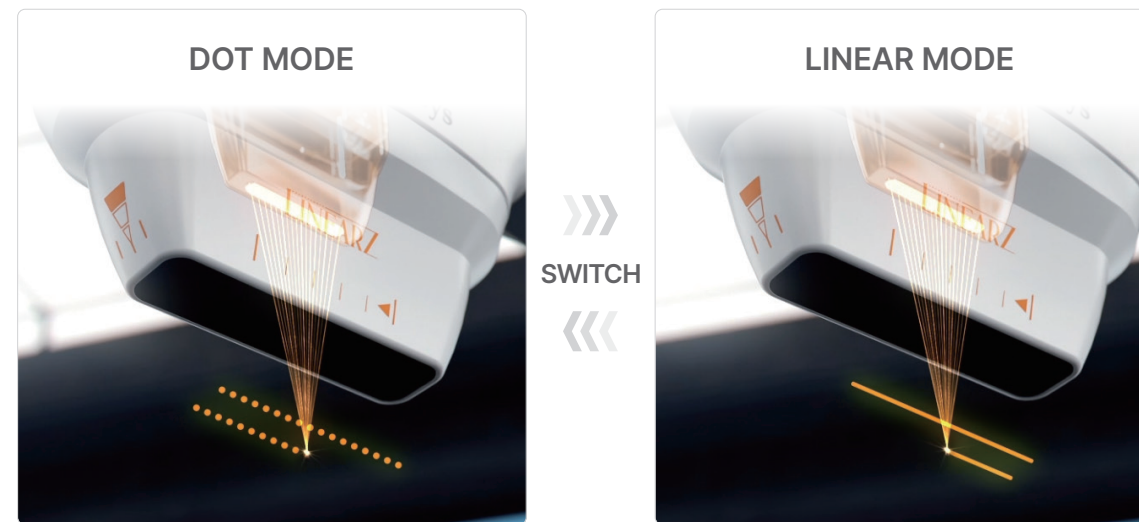

Energy is delivered in a dot pattern with 0.1mm spacing increments.

Energy is delivered seamlessly without spacing.

DWR System (Degassed Water Replacement)

Circulates internally stored degassed water to prevent cartridge overheating, enabling stable and continuous treatments through rapid reactivation.



Details Made For You

Handpiece designed for safe and comfortable treatments.



User-Friendly UI

Allows real-time monitoring of treatment information such as Mode & Depth, shot count, cartridge temperature, and transducer position.